


URL Calls and Sharedlinks

Shared links (hanging protocols)

Suppose you have arranged a case in the viewer and want to send it to partner, you can save the state and generate a link to access the state via

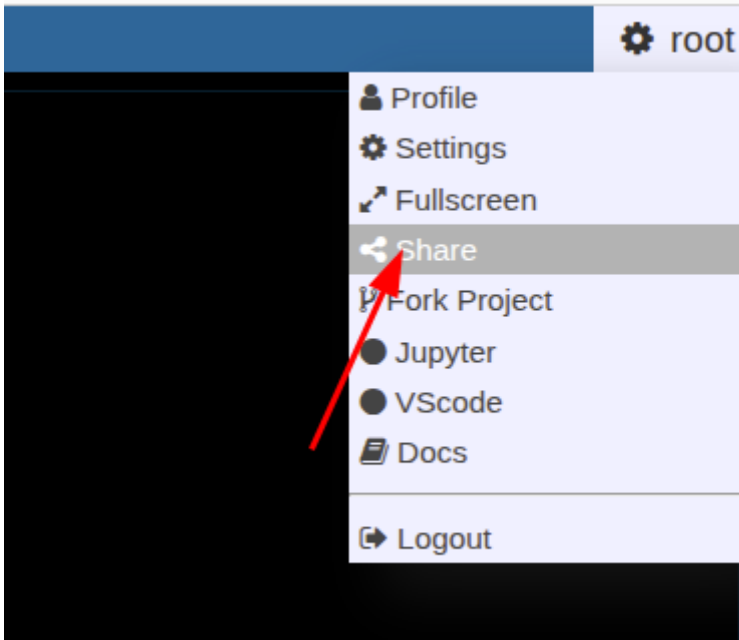

If you have called a link, you can modify and use again the "share" button to overwrite the old link (or create a new one).

URL calls

Suppose you have a NORA project with configured autoloader and want to call a specific case from an external application. You can call NORA via

```
https://nora.ukl.uni-  
freiburg.de/godzilla/index.php?asuser=root&project=WHATEVER&call={"pid":"yourpatientid","sid":"your  
studyid","selmode":"yourmode","preset":"yourpresetname","user":"username","token":"yourtoken"}
```

The only mandatory options are "user" and "token" for authentication (if you are not already authenticated). Get the token from your profile. Note that you need to create your own password (LDAP authentication is not possible).

passwd:

click to change

token:

click to copy to clipboard

Styling theme:

classic blue ▼

Size of local hard disk cache, use this if you are remote (in MB):

1

Revision #5

Created 5 June 2025 07:38:35 by reisertm

Updated 6 June 2025 10:00:48 by reisertm